

# A3-AI intelligent skin tester

LEAD TECHNOLOGY CREATE BEAUTY

# A3 AI intelligent skin tester

A3 AI intelligent skin tester: With a research and development focus on solving skin problems, it integrates eight spectral imaging technologies and can professionally and objectively analyze seventeen problems of facial skin with flexible operations. The original intention of the research and development is to take photos and analyze reports with just one click, making it more convenient to operate.



## Supports 19 languages



Traditional Chinese



English



French



German



Japanese



Korean



Spanish



Portuguese



Italian



Russian



Dansk



Dutch



Polish



Turkish



Arabic



Tiếng Việt



Indonesian



Thai

【A3 AI intelligent skin tester】

Lead technology create beauty

## Adapt to the scene





【A3 AI intelligent skin tester】

# Catalogue

01

FUNCTION  
DEMO

02

PARAMETERS

03

EIGHT  
SPECTRAL  
IMAGE  
ANALYSIS

04

22  
DETECTION  
FUNCTION

05

MICROSCOPIC  
DETECTION  
SKIN DETAILS

06

OUR  
SERVICE

【A3 AI intelligent skin tester】


Lead technology create beauty



A3 AI intelligent skin tester

01

FUNCTION  
DEMO

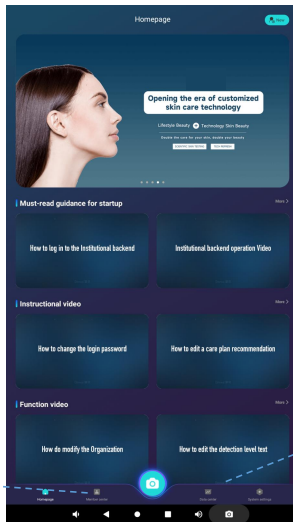

Banner

Must-read guidance  
for startup

Instructional video

Function video

Member Center

Data Center


Click on the setting to adjust  
the parametersters

【A3 AI intelligent skin tester】

Lead technology create beauty

Analysis of 4 major symptoms

30+ detection dimensions

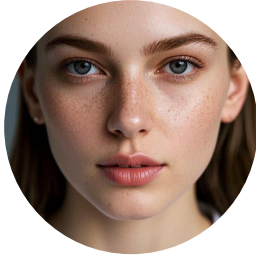

### Analysis of aging



Forehead lines



Dorsal nasal lines



Lines around the eyes



Crow's feet



Nasolabial folds

### Sensitive analysis



Acne



Redness



Acne rosacea



Barrier

### Pigment analysis



Mole



Freckles



Acne marks



Spots

### Skin quality analysis



Pores



Porphyrin



Wrinkle



Moisture

【A3 AI intelligent skin tester】

Lead technology create beauty

## Function demo

THREE DIMENSIONAL SHOOTING



TWO SKIN TESTING MODES



SEVENTEEN TESTING INDICATORS



THREE COMPARATIVE ANALYSIS MODES



THREE REPORTING MODES



A3 AI intelligent skin tester

EIGHT SPECTRAL IMAGING



CLOUD DATA STORAGE



CUSTOMER FILE MANAGEMENT



ADJUSTING SKIN TEST RESULTS




BACKEND SYSTEM MANAGEMENT




【A3 AI intelligent skin tester】

Lead technology create beauty


## Three dimensional shooting

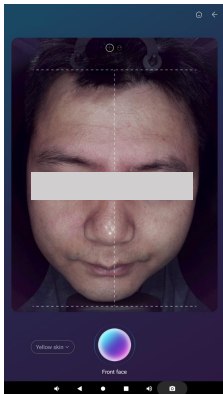

Front face

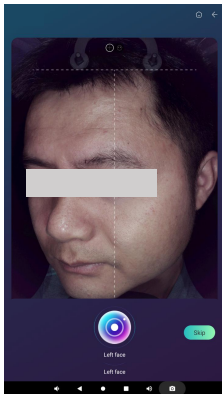


Left face

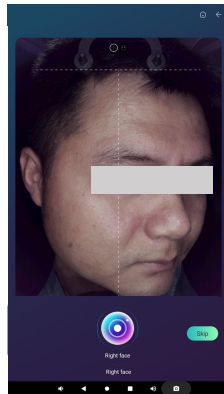

Right face

## Four skin tones available

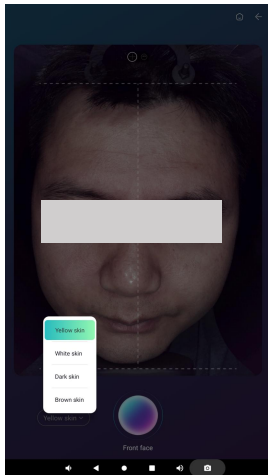


**Yellow Skin**

**White skin**

**Dark skin**

**Brown skin**

- Eight spectral images

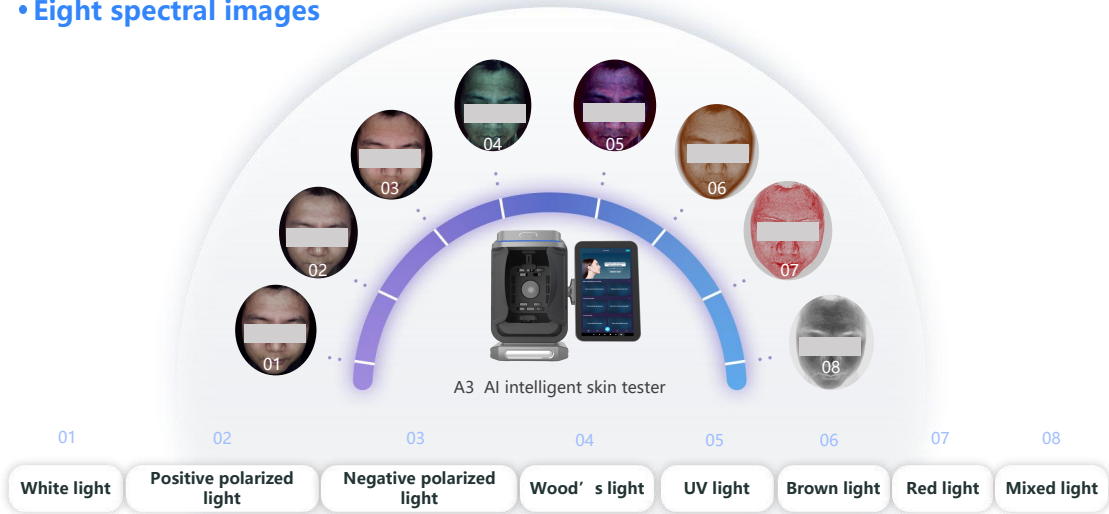






## Functional coverage

- Microscopic imaging of skin details presentation



**Handheld skin microimager**

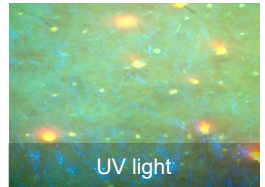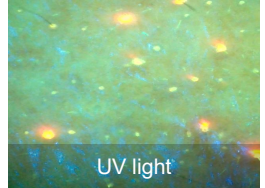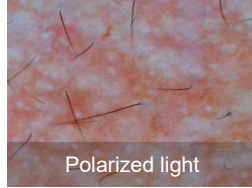

**Multiple light source targeting analysis**

## Seventeen testing indicators



Hydration



Pores



Blackhead



Sebum



Skin Barrier



Acne



Wrinkle



Mixed spot



Superficial pigment



A3 AI intelligent skin tester



Porphyrin



Collagen



Fluorescent agent



Deep pigment



Brown pigment



Heat Map of Sensitivity



Heat Map of Pigment




Red Map of Sensitivity

## Three comparison modes




Horizontal comparison

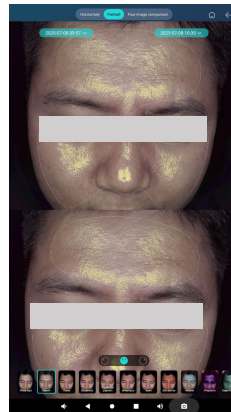

Vertical comparison


## Three comparison modes

Comparison of the effects of single indicators before and after.

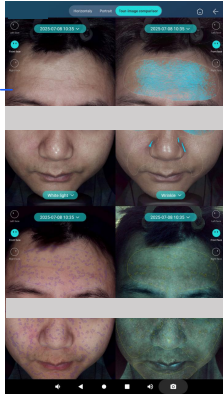


Comparison of four figures

Analysis and comparison of multiple problem indicators before and after skin care.

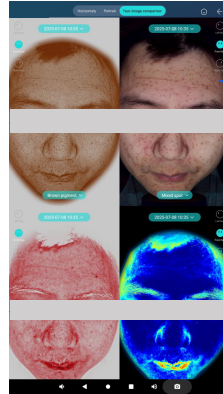

Comparison of four figures

## Three reporting modes-single independent report 1

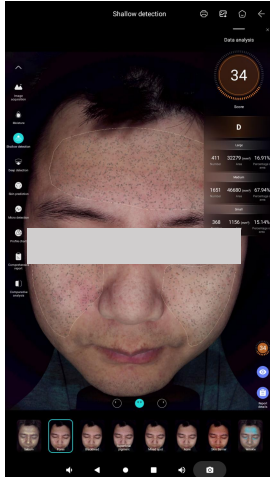

### Introduction to data analysis

- 1.Score
- 2.Level
- 3.Number
- 4.Area
- 5.Percentage of area

Let consumers accurately understand their skin problems

Accurately quantify the underlying effect and empower doctors to treat


### Data score



It is divided into 5 levels according to the skin condition from high to low and marked with different colours.



A Green B Blue C Yellow D Orange E Red

### Multidimensional perspective

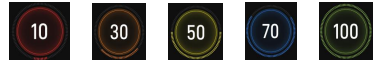

Multi-angle comprehensive image display

## Three reporting modes-single independent report 1

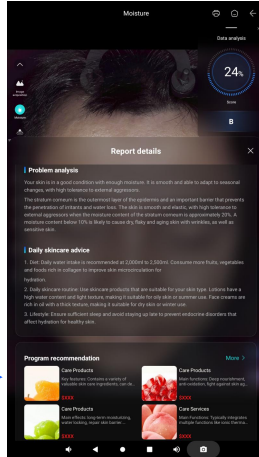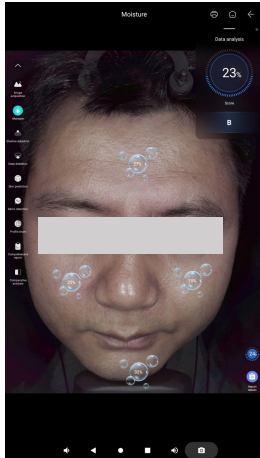


Problem analysis

Nursing advice

Scheme recommendation



## The classification of 16 skin types in Baumann Skin Typing System

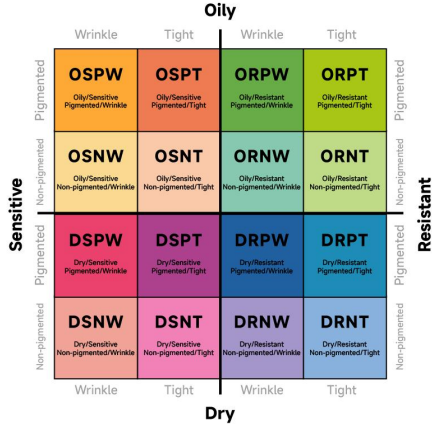

16Type	Explanation	16Type	Explanation
DSPT	Dry, Sensitive, Pigmented, Tight	OSPT	Oily, Sensitive, Pigmented, Tight
DSPW	Dry, Sensitive, Pigmented, Wrinkle	OSPW	Oily, Sensitive, Pigmented, Wrinkle
DSNT	Dry, Sensitive, Non - Pigmented, Tight	OSNT	Oily, Sensitive, Non - Pigmented, Tight
DSNW	Dry, Sensitive, Non - Pigmented, Wrinkle	OSNW	Oily, Sensitive, Non - Pigmented, Wrinkle
DRPT	Dry, Resistant, Pigmented, Tight	ORPT	Oily, Resistant, Pigmented, Tight
DRPW	Dry, Resistant, Pigmented, Wrinkle	ORPW	Oily, Resistant, Pigmented, Wrinkle
DRNT	Dry, Resistant, Non - Pigmented, Tight	ORNT	Oily, Resistant, Non - Pigmented, Tight
DRNW	Dry, Resistant, Non - Pigmented, Wrinkle	ORNW	Oily, Resistant, Non - Pigmented, Wrinkle

## Three reporting modes-Comprehensive analysis report 2



Personal Information

16-type skin classification  
and comprehensive score

On the left is the skin  
water content and on  
the right is a single  
report.

Single indicators below C will be  
realised in the form of radar  
charts.

Comparative analysis,  
only compare and  
analyse with the results  
of the last shooting.

Shallow test

Deep test




## Three reporting modes-comprehensive analysis report 2



### Overall recommendation

Real-time analysis of problems and nursing suggestions based on Bauman's 16 types of skin classification


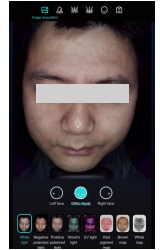


H5 Mobile phone report

### Download the report

Scan the code on your mobile phone to get the report

## Three kinds of reporting modes Innovation Independent Editing Report 3

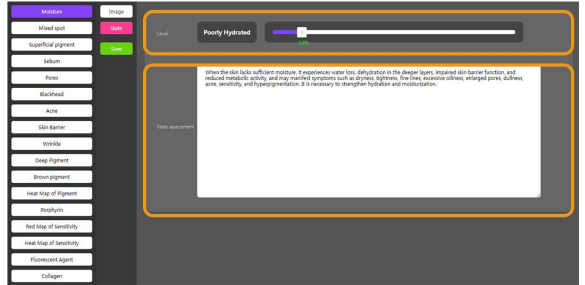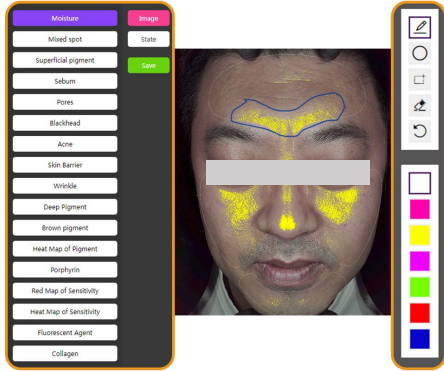


The background automatically selects the image to be edited for annotation

Customize any test results you want

You can debug results for each indicator

Free copy editing


## Skin prediction

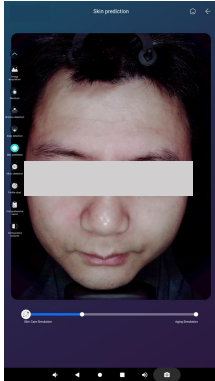

Nursing simulation


Deeply predict the future of the skin and awaken customers' desire for young skin.

According to the current skin condition of the customer, through the training of AI large model, simulate the skin condition after customer care and the aging situation of different age groups.

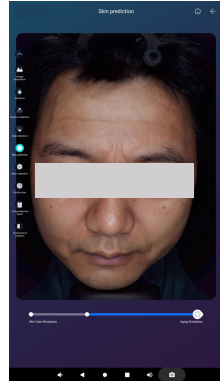


Aging simulation

## Parameter adjustment

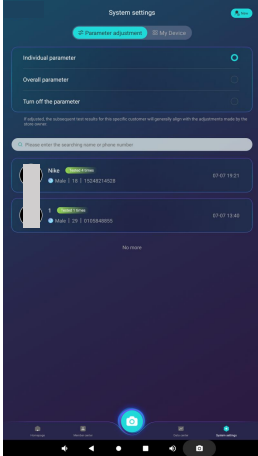


Parameter  
adjustment

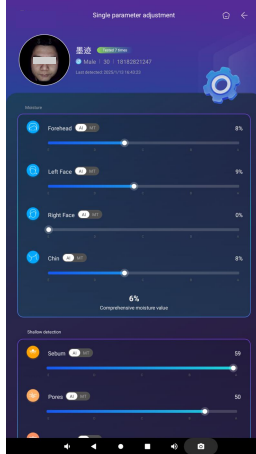


Data result  
optimization  
manual debugging

It can be adjusted as a whole

## My device

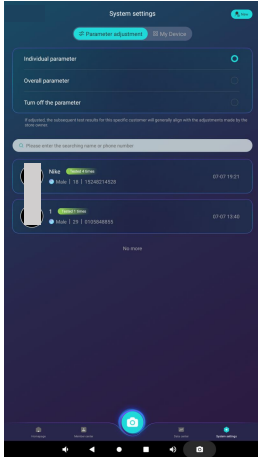




The language  
can be adjusted

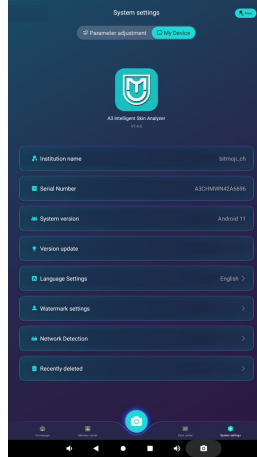

Serial number

You can find the instrument  
problem through the  
background and solve it.

## Watermark function

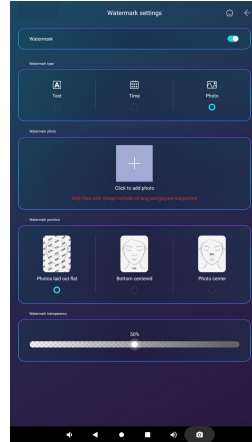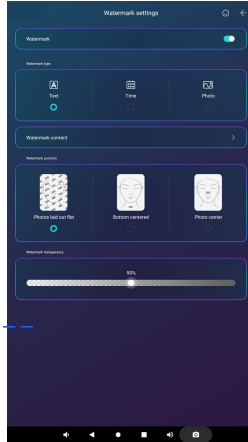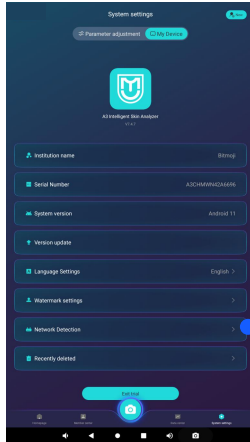

Customizable  
text watermark,  
image  
watermark

## Birthday reminder function

Create New Profile

Username Please enter user name


Gender ☐ Male ☐ Female

Birthday Select date of birth >

Phone number Please enter the phone number

Consultant Please enter the consultant's name

Save

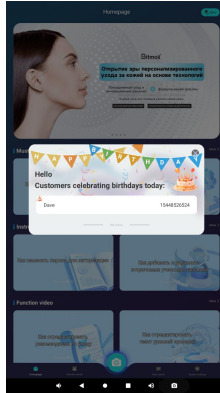

Automatically send out the list of customers whose birthdays are on that day, in order to enhance customer retention.

## Image export function



### Exporting images

#### Hint

Save to the following path:  
/storage/emulated/0/A3Export/Pictures/116207962/20250716\_202231\_White\_Light\_2.jpg

Confirm

Turn on

This figure is exported.

#### Select photo to be exported

☐ Select all

Image



White\_Light



Negative\_photo



Positive\_photo



White\_Light



UV\_Light



Red\_Light



Brown\_photo



Annotated photo

☐ Select all



Baseline



Acne



Collagen



Deep\_pigment



Fluorescent\_Agent



Heat\_Map\_of\_Deep\_pigment



Heat\_Map\_of\_Pigment



Sebum



Wrinkle



Pore



Moisture



Superficial\_pigment




Export

Export multiple images




## Member Center

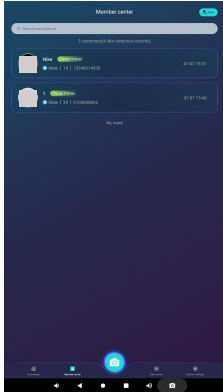

File Management

The test items with scores below 60 will be listed here.

Bauman 16-type classification

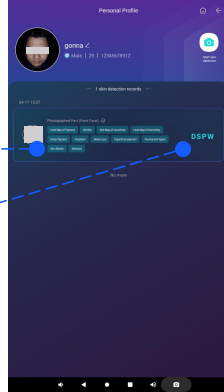

Profile

Cloud storage

One click search

Number of skin tests

## Data Center

### Choose the year and the month

The number of people entering the store in a month or a year can be displayed more intuitively through the line chart

### Ratio of men to women

He men and women who enter the store can see the ratio of male and female women who enter the store through the data form more clearly and intuitively

### Problem proportion

In the form of data, the current population is divided into five levels with different colors to represent the data in each indicator



### Available year

### Age distribution

Through the age division of the number of people entering the store in the form of data.

## Backend management



Recording  
nursing  
project freely

Unified  
management  
of multiple  
equipments

Customer  
profile  
management  
in real-time

Review  
and edit  
detection  
record in  
real-time

【A3 AI intelligent skin tester】

Lead technology create beauty



**A3 AI intelligent skin tester**

02

**HARDWARE  
PARAMETERS**

## Parameters

Device Name	A3-AI intelligent skin tester
Model	ZMLH-A3
Pixel	Industrial grade 36 Million pixel
Spectrum	RGB white light; Positive polarized light; Negative polarized light; Wood's light; UV light;Red light; Brown light; Mixed light
Shading Method	Semi-open
Product Material	Industrial Grade ABS
Power Supply	AC 100-260V, 50/60Hz
Product Dimensions	Unfolded: 600 x 415 x 504mm; Folded: 316 x 415 x 504mm
Central Processor	Rockchip RK3568 ; Quad-core Cortex-A55
Motherboard	R10-S6810 Pmotherboard,running Android 11
Main Memory	Dual-channel LP DD R4, 4G
Storage	MMC 5.1, 32G

## Parameters

Device Name	A3-AI intelligent skin tester
Operating Method	multi-point capacitive touch
Screen Features	Foldable 90 degrees
WIFI	Built-in dual-band WIFI (2.4G, 5G)
HDMI	1 port
USB	2 ports
Screen Ratio	16:9
Screen Size	15.6 inches
Screen Resolution	1920*1080
Net Weight	11.8 KG
Gross Weight	16.6 KG
Packaging Dimensions	595 x 570 x 620mm
Packaging Material	Corrugated paper + pearl cotton
Other Accessories	Handheld skin analyzer, power cord

## Hardware parameters



【A3 AI intelligent skin tester】

Lead technology create beauty



A3 AI intelligent skin tester

03

EIGHT  
SPECTRAL  
IMAGE  
ANALYSIS



## Eight spectral image analysis



White light


### THEORY

Visible spots and other blemishes on the skin surface (acne, spots, wrinkles, pores, etc.) under natural light sources, which are mainly used as the basis for other spectral image comparison.


Acne



Dark spots

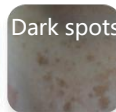


Wrinkle

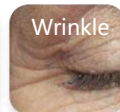

Pore



## Eight spectral image analysis

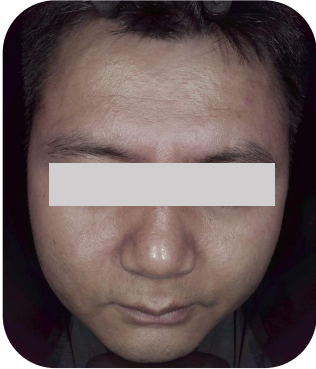

Positive polarized light

### THEORY

Positive polarized light can improve the clarity of superficial texture, magnify local details, so as to clearly observe the smoothness of skin, fine lines and wrinkles and bumps (wrinkles, pores, Acne scars, Acne, etc).



## Eight spectral image analysis




Negative polarized light

### THEORY

Using negative polarized technology to filter out the refracted light on the skin surface, so that you can clearly examine the light brown, tan, dark brown, light yellow or dark red skin lesions; It can distinguish the condition of capillaries, facial acne, uniformity skin and other skin problems.

Dark spots

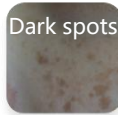


Pigment



Acne



## Eight spectral image analysis

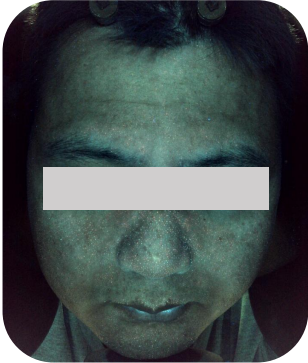

Wood' s light

### THEORY

Wood' s light can detect deep pigments in dermis. The principle behind this is that melanin does not fluoresce after exposure to ultraviolet radiation, allowing melanin to stand out more clearly with stronger contrast.

Spots

Spots

Fluoresce-  
nce

## Eight spectral image analysis



UV light

### THEORY

Under UV light source, the content and distribution of the purple pigment bilirubin are displayed clearly through fluorescence, which can be used for the auxiliary diagnosis and efficacy observation of pigmentary dermatoses, pore issues, skin infections, and porphyria.



Porphyrin



Fluorescence

## Eight spectral image analysis



Brown light

### THEORY

The position, area, shape, and severity of subcutaneous facial UV spots are processed by using RBX light source technology, which demonstrate skin damage from UV radiation and the accumulation of subcutaneous melanin.



## Eight spectral image analysis




Red light

### THEORY

Using the RBX light source technology, it is applied for analyzing the subcutaneous hemoglobin on the face as well as inflammatory pigmentation (such as sensitive skin, lesions, acne, red blood vessels, etc, which are skin inflammatory conditions)

Acne

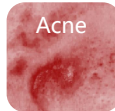


Blood  
streak



Mole



## Eight spectral image analysis

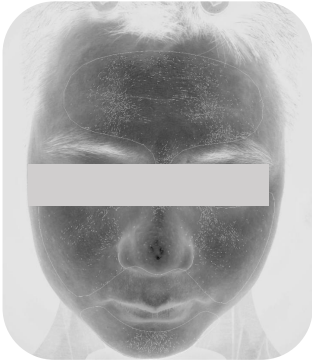

Mixed light

### THEORY

Skin texture roughness and collagen loss were revealed by polarizing analysis.

Rough  
texture

Wrinkle



【A3 AI intelligent skin tester】

Lead technology create beauty




**A3 AI intelligent skin tester**

04

22

**DETECTION  
FUNCTION**

## 22 Detection function-Moisture test report 1

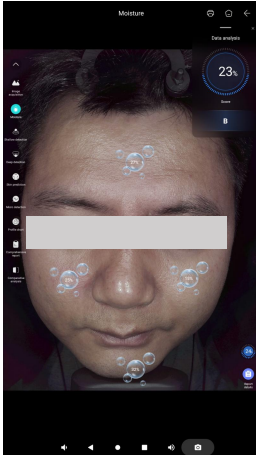

### Zoned moisture detection

Partition detection of skin and facial moisture

### Skin water content score



The moisture content of the skin is classified into five levels from low to high, and each level is marked with a different color.



E



D



C




B



A

## 22 Detection function-Moisture test report 2

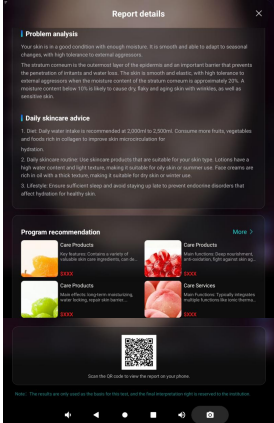


### Problem analysis


**Nursing advice**

### Scheme recommendation


Scan the code to view the report

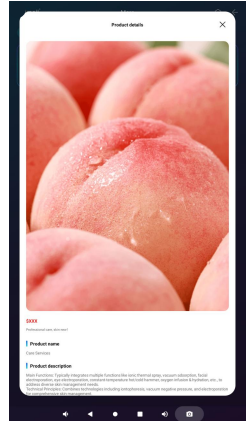

## 22 Detection function-Under - skin perspective

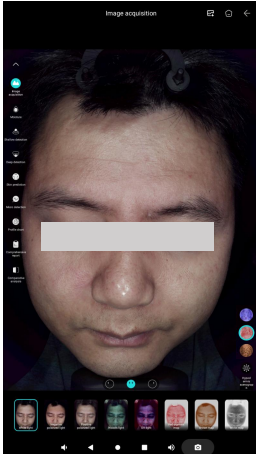


Under - skin  
perspective

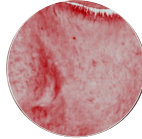

Clicking the "Muscle Layer Perspective" button on the right side will reveal three icons: the UV light map, the melanin map, and the brown map. You can freely move the perspective lens to compare issues such as pigmentation, red blood vessels, and skin inflammation on the superficial layer and the base layer of the skin.

## 22 Detection function


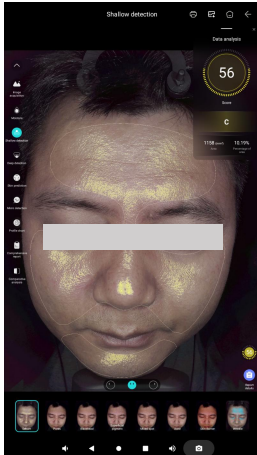


### Sebum

#### Image Analysis

- ❑ The oil secretion of the skin surface can be checked under positive polarized light source.
- ❑ The algorithm displays areas of the skin with active oil secretion through yellow ,through the form of data, you can see the oiliness of facial skin more clearly and intuitively.
- ❑ Excess oil is one of the factors that trigger acne growth, so please take good oil control care if you have acne.

## 22 Detection function

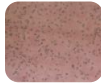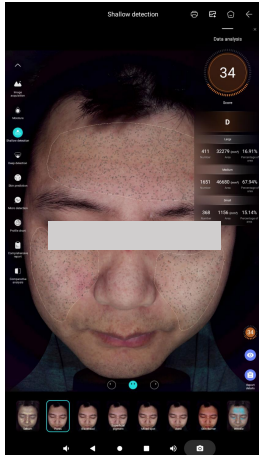


### Pores

#### Image Analysis

- ❑ Under negative polarized light source, it is possible to check if enlarged pores have formed on the skin surface.
- ❑ The algorithm uses the RBX technology to process and present the areas with large pores in the skin as dark gray small circles; it can also show the situation of already enlarged pores in the facial skin in a clearer and more intuitive data form.
- ❑ Pore clogging refers to the pores on the surface of the skin being blocked, which prevents sebum from being discharged normally, accompanied by the accumulation of stratum corneum and dirt. this phenomenon usually manifests itself in the form of blackheads, whiteheads or acne, and in severe cases may lead to skin problems such as acne and folliculitis.

## 22 Detection function



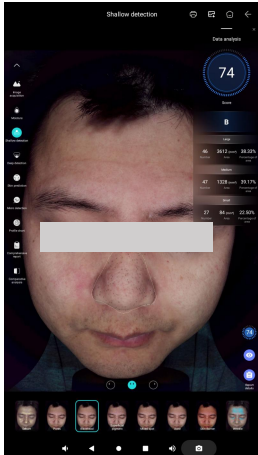


### Blackhead

#### Image Analysis

- ❑ Under negative polarized light source, you can see blackheads formed by pores clogged by oil in the T-zone.
- ❑ The algorithm uses the RBX technology to process and present the area of blackheads in the T-zone as small dark gray circles; it enables a clearer and more intuitive view of the blackheads on the nose area through the data.
- ❑ Blackheads are formed by excess oil accumulation in the nose area of the skin and air oxidation. Areas with large pores are more likely to accumulate and store oil and dust in the air, so it is necessary to clean and moisturize in time to reduce the formation of large pores.

## 22 Detection function

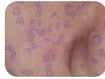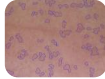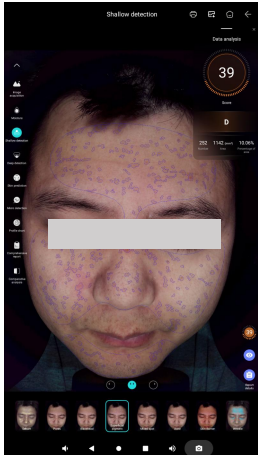




### Superficial Pigment

#### Image Analysis

- ❑ Superficial pigment refers to pigmentation that has formed on the superficial layer of the skin, including: acne scars, spots, inflammatory pigmentation, etc.
- ❑ The coverage of pigmentation may exist in both deep and shallow layers. you can compare the image with the deep pigment image. If the shallow layer shows pigmentation but the deep layer shows no pigmentation, it means that the pigment is only deposited in the superficial layer of the skin.
- ❑ The algorithm marks the pigmented area with a purple polygon curve, and the shallow pigment can be seen more clearly and intuitively through the form of data.



## 22 Detection function

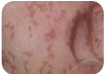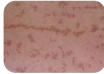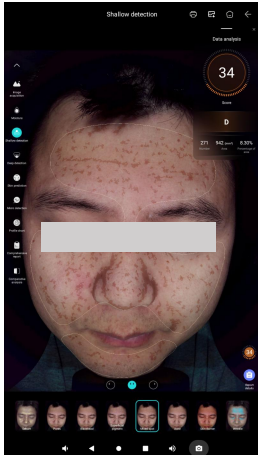


### Mixed Spot

#### Image Analysis

- ❑ Under negative polarized light, we can see the distribution of mixed spots on the skin surface. The algorithm identifies the facial complex spot area and marks it with a brown block.
- ❑ The mixed spot map shows skin pigmentation such as melasma, age spots, and freckles. Melasma is a darker patch on the skin that can appear brown, black, or dark brown. Melasma may expand over time, especially if daily sun protection and skin care are not taken care of. Some melasma may be slightly raised and feel slightly convex to the touch.
- ❑ The algorithm marks the mixed spot area with brown color blocks, and the mixed spot situation can be seen more clearly and intuitively through the form of data.

## 22 Detection function

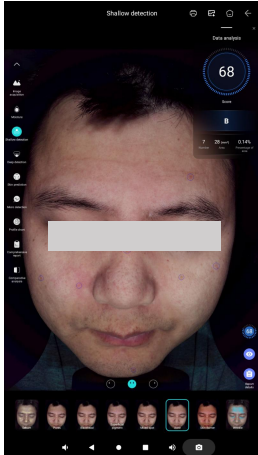


### Acne

#### Image Analysis

- ❑ Look at the distribution of skin acne and superficial redness under negative polarized light.
- ❑ When pores are clogged with oil and dust, it is easy to fester or form inflammation, which will then turn into acne and acne.
- ❑ The algorithm identifies the distribution area of facial acne and marks it with blue circles. The more the number and the more obvious the redness of the skin, the more serious the skin acne problem is, and the skin needs to be oil-controlled to unclog the pores and eliminate inflammation. You can see the acne situation more clearly and intuitively through the form of data.

## 22 Detection function


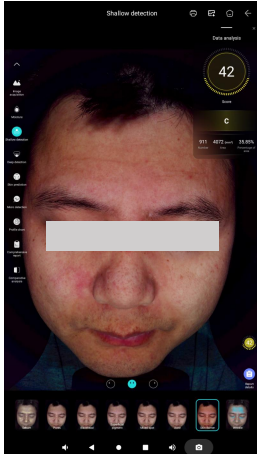


### Skin Barrier

#### Image Analysis

- ❑ We can check the skin barrier health under negative polarized light source.
- ❑ The barrier image shows the skin redness problem and the distribution of red blood streaks. The formation of red blood streaks is mainly due to the damage of keratin, the weakness of the epidermis, and the long-term damage of the capillary position, which leads to vascular dilation and congestion.
- ❑ The redder the area, the more severe the damage to the skin barrier. This can be used as a reference for judging the sensitivity of the skin and the presence of inflammatory areas. The damaged condition of the barrier can be seen more clearly and intuitively in the form of data.

## 22 Detection function

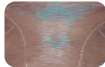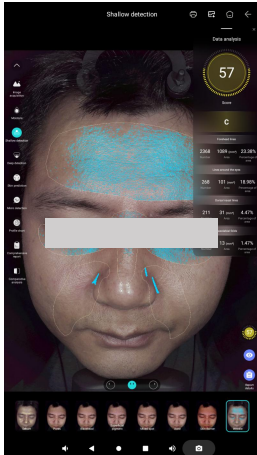


### Wrinkle

#### Image Analysis

- ❑ The texture of the skin surface can be viewed under a positive polarized light source.
- ❑ The wrinkle image shows the roughness of the skin texture, such as large pores, dry lines, fine lines, and wrinkles. It can be used as a reference for judging the fineness of the skin and the loss of collagen.
- ❑ The algorithm identifies the patterns on the facial skin and marks the distribution of skin wrinkles in five areas (forehead wrinkles, nasal bridge wrinkles, eye area wrinkles, crow's feet wrinkles, and laugh lines) with Lake blue dotted lines. The more discontinuous the lines, the rougher the skin. This way, the wrinkles can be seen more clearly and intuitively in a data form.

## 22 Detection function



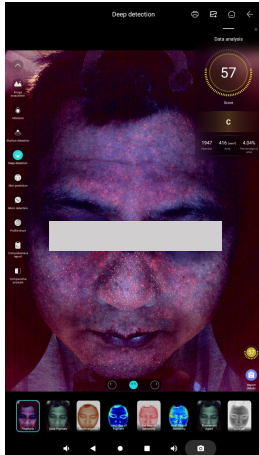


### Porphyrin

#### Image Analysis

- ❑ The brick-red fluorescent spots in the picture are propionibacterium acnes and malassezia. these two bacteria will aggravate the occurrence of skin acne, so they can be used as a basis for judging skin acne. Through the form of data, the situation of porin can be seen more clearly and intuitively.
- ❑ The living environment of propionibacterium acnes and malassezia must have oil, so they can be used as a basis for judging the accumulation of oil in skin pores.

## 22 Detection function

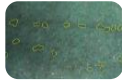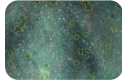
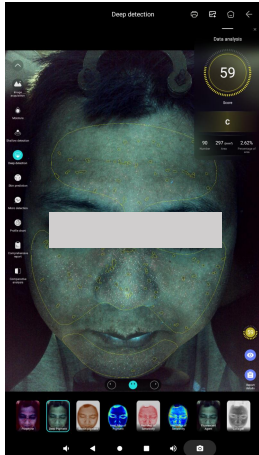


### Deep Pigment

#### Image Analysis

- ❑ The bottom part of the figure shows wood's light. The yellow area represents the comprehensive color spot area identified by the algorithm and is marked with polygonal curves. This way, the situation of deep pigmentation can be seen more clearly and intuitively in the form of data.
- ❑ The dark (black, brown) block or dot skin that appears on the face is a display of skin pigmentation (such as melasma, freckles, malar spots, inflammatory pigmentation, acne marks, hemoglobin aggregation, etc.).
- ❑ The pigmentation in the deep layer of the skin can be compared with the sensitivity to determine whether it is an inflammatory hemoglobin accumulation or a spot problem.

## 22 Detection function



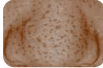
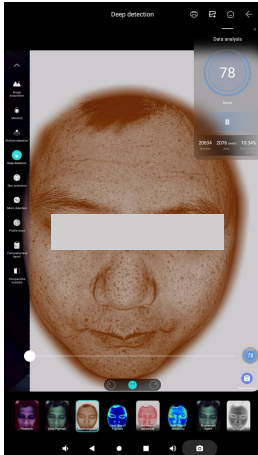


### Brown Pigment

#### Image Analysis

- ❑ The depth of the overall brown color of the skin is mainly related to the skin color. People with darker skin or more hemoglobin have darker overall pigmentation.
- ❑ The areas with heavier pigmentation in the image are mostly those with higher pigment concentration density.
- ❑ Through the form of data, the brown pigment can be seen more clearly and intuitively.

## 22 Detection function

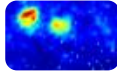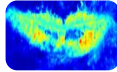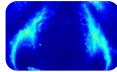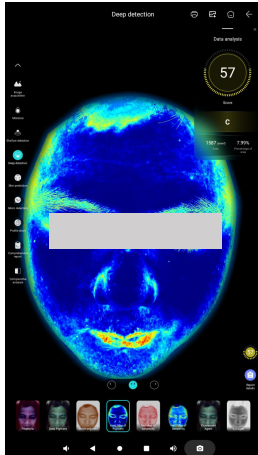




### Heat Map of Pigment

As shown in these images

- ❑ Heat map of pigment can check the distribution of pigment deep in the skin.
- ❑ The algorithm identifies the distribution of pigments on the face and presents it in the form of a heat map. Different colors are used to represent the distribution of spots, moles, and scars visible to the naked eye under negative polarized light. Red indicates severe skin pigmentation, yellow for medium, green for lighter skin, and blue for normal skin. The situation of pigmentation can be seen more clearly and intuitively through the form of data.
- ❑ Pigment production mechanism: the body's own regulation, physical or chemical factors stimulate melanocytes, so that their number increases and activity increases. The melanin produced cannot be completely removed with the stratum corneum and blood circulation, and finally deposits in the local skin.



## 22 Detection function


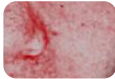
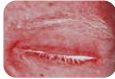
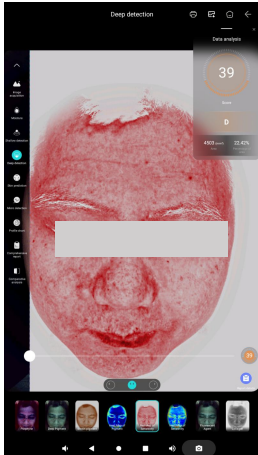


### Red Map of Sensitivity

#### As shown in these images

- ❑ In the red map of sensitivity, we can observe the redness of the skin's surface layer and the distribution of red blood vessels.
- ❑ The picture clearly shows the distribution of red blood vessels, indicating that the skin is thin and sensitive, and requires proper protection and care.
- ❑ The shade of the background color in the red photo is related to the overall skin tone. Those with less hemoglobin will have a lighter color.
- ❑ The areas with a higher concentration of red color indicate a more concentrated accumulation of hemoglobin in the skin, which can be used as a reference for judging the sensitivity of the skin and the presence of inflammatory areas. The degree of skin sensitivity can be clearly and intuitively seen through the data form.

## 22 Detection function

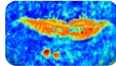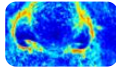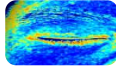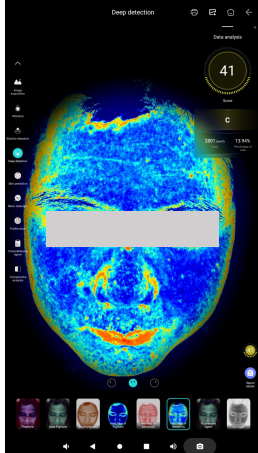


### Heat Map of Sensitivity

#### As shown in these images

- ❑ The heat map of sensitivity represents skin sensitivity. When the skin shows significant redness and thinning of the stratum corneum, it becomes more susceptible to external stimuli and damage, leading to issues such as dryness, sensitivity, and redness.
- ❑ The heat map of sensitivity is based on the distribution of subcutaneous capillaries, with areas of greater sensitivity having more capillaries. Visible redness and acne on negative polarized light images indicate areas of severe sensitivity.
- ❑ The algorithm uses different colors to indicate varying degrees of sensitivity and their distribution on the skin. Areas with severe sensitivity are shown in deep red, including the lips; medium sensitivity is represented in yellow, mild sensitivity in green, and normal skin appears in blue. The sensitivity is more clearly and intuitively reflected in the form of data.

## 22 Detection function


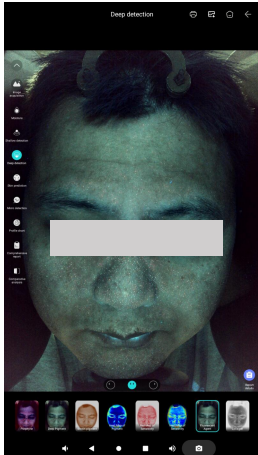


### Fluorescent Agent

#### As shown in these images

- ❑ Fluorescent agent and pigments might both appear in facial imaging. To assess the fluorescent agent, focus specifically on the fluorescence response.
- ❑ The difference between fluorescent agent and porphyrins is as follows porphyrins exhibit brick-red fluorescent agent, while fluorescent agent display intense blue light and usually appear as large, sheet-like areas.
- ❑ The difference between fluorescent agent and facial .
- ❑ Dust is as follows facial dust appears as white, bright, floating, and short, wispy lines on the surface, while fluorescent agent typically display bright colors and are often more diffuse or spread over larger areas.

## 22 Detection function

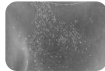
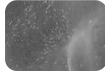


### Collagen

#### As shown in these images

- ❑ Under the mixed light image, we can observe the loss of collagen on the skin surface.
- ❑ The mixed light image shows a situation where the skin texture is rough. For example: large pores on the skin, fine lines and wrinkles. It can be used as a reference for judging the skin's smoothness and the loss of collagen.
- ❑ The white dotted lines in the mixed light image indicate the loss of skin collagen and the rupture of elastic fibers. This way, the loss of collagen can be clearly and intuitively observed in a data form.

## Profile chart

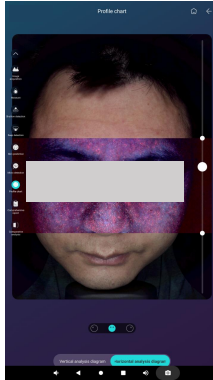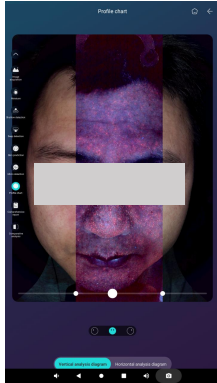

### Profile chart

Through white light, negative polarized light and UV light source comparison, multi-dimensional, deep analysis of skin problems.

【A3 AI intelligent skin tester】

Lead technology create beauty




**A3 AI intelligent skin tester**

05

**MICROSCOPIC  
DETECTION  
SKIN DETAILS**

## Microscopic detection - Skin Details

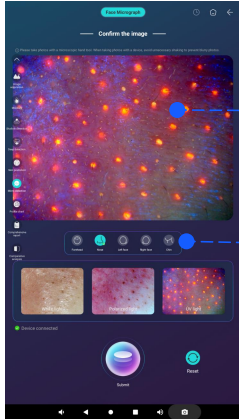

Local microscopic display



Location selection, partition detection.



## Microscopic detection - Skin Details



Independent  
record tracking  
for each  
customer

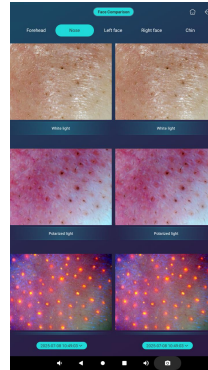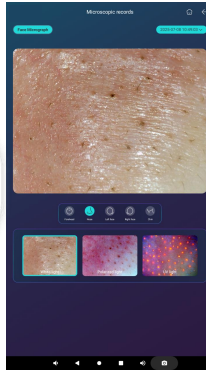

Comparison  
display for the  
same area



【A3 AI intelligent skin tester】

Lead technology create beauty




**A3 AI intelligent skin tester**

06


OUR  
SERVICE

## Our services

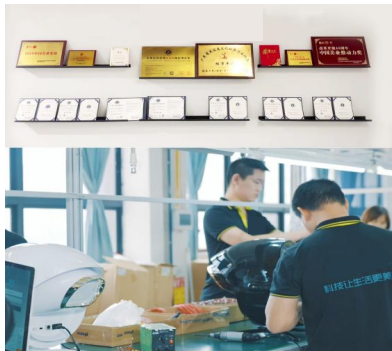

Meticulous  
craftsmanship

Original  
manufacturer

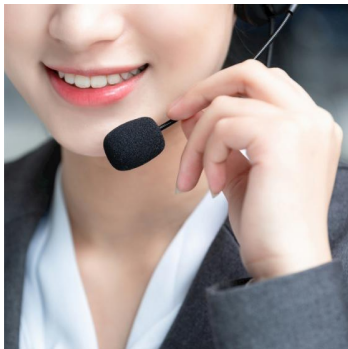

Specialist  
service

After-sales  
service



Training and  
guidance

Specialized  
teaching

# Our services

【A3 AI intelligent skin tester】

Educational materials, user manuals, and instructional videos are all provided.



Welcome to the world of Moji!  
Explore skin secrets with Moji!

AI Intelligent skin tester  
User manual

AI intelligent skin tester (2MLH-A3) product list

Model	Color	Quantity	Accessories	Remarks
2MLH-A3	White	1	1	
2MLH-A3	White	1	1	
2MLH-A3	White	1	1	
2MLH-A3	White	1	1	
2MLH-A3	White	1	1	
2MLH-A3	White	1	1	

Model 2MLH-A3

QR code

AI intelligent skin tester MAX version |  
operation explanation video



AI intelligent skin tester Plus version |  
operation explanation video



AI intelligent skin tester Plus version |  
Unboxing and installation video



Hyperbubble skin rejuvenator (Unboxing video)



Hyperbubble skin rejuvenator (Function introduction)



Hyperbubble skin rejuvenator (Method explained)



Hyperbubble skin rejuvenator (How to operate it)



---

【A3 AI intelligent skin tester】

Lead technology create beauty

# Welcome you join us!

THANK YOU FOR WATCHING