

B6 AI Intelligent skin analyzer

B6 AI Intelligent skin analyzer : With a research and development focus on solving skin problems, it integrates eight spectral imaging technologies and can professionally and objectively analyze seventeen problems of facial skin with flexible operations. The original intention of the research and development is to take photos and analyze reports with just one click, making it more convenient to operate.



B6 AI Intelligent skin analyzer

Supports 19 languages



Traditional Chinese



English



French



German



Japanese



Korean



Spanish



Portuguese



Italian



Russian



Dansk



Dutch



Polish



Turkish



Arabic



Tiếng Việt



Indonesian




Thai

【B6 AI Intelligent skin analyzer】

Lead technology create beauty

Adapt to the scene

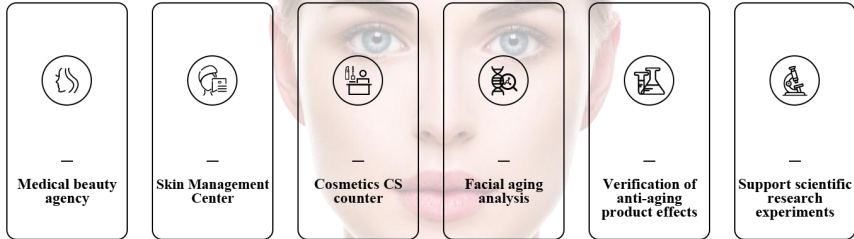

【B6 AI Intelligent skin analyzer 】

Catalogue

01

**FUNCTION
DEMO**

02

PARAMETERS

03

**EIGHT
SPECTRAL
IMAGE
ANALYSIS**

04

**24
DETECTION
FUNCTION**

05

**MICROSCOPIC
DETECTION
SKIN DETAILS**

【B6 AI Intelligent skin analyzer】

Lead technology create beauty




B6 AI Intelligent skin analyzer

01

FUNCTION
DEMO

【B6 AI Intelligent skin analyzer】

Lead technology create beauty

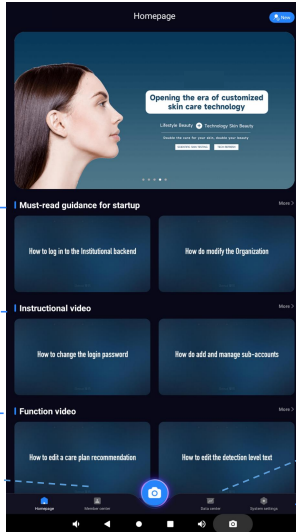

Banner

Must-read
guidance for startup

Instructional video

Function video

Member Center

Data Center


Click on the setting to adjust
the parametersters

【B6 AI Intelligent skin analyzer】

Lead technology create beauty

Analysis of 4 major symptoms

30+ detection dimensions

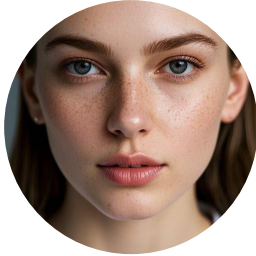

Analysis of aging



Forehead lines



Dorsal nasal lines



Lines around the eyes



Crow's feet



Nasolabial folds

Sensitive analysis



Acne



Redness



Acne rosacea



Barrier

Pigment analysis



Mole



Freckles



Acne marks



Spots

Skin quality analysis



Pores



Porphyrin




Wrinkle




Moisture


Function demo

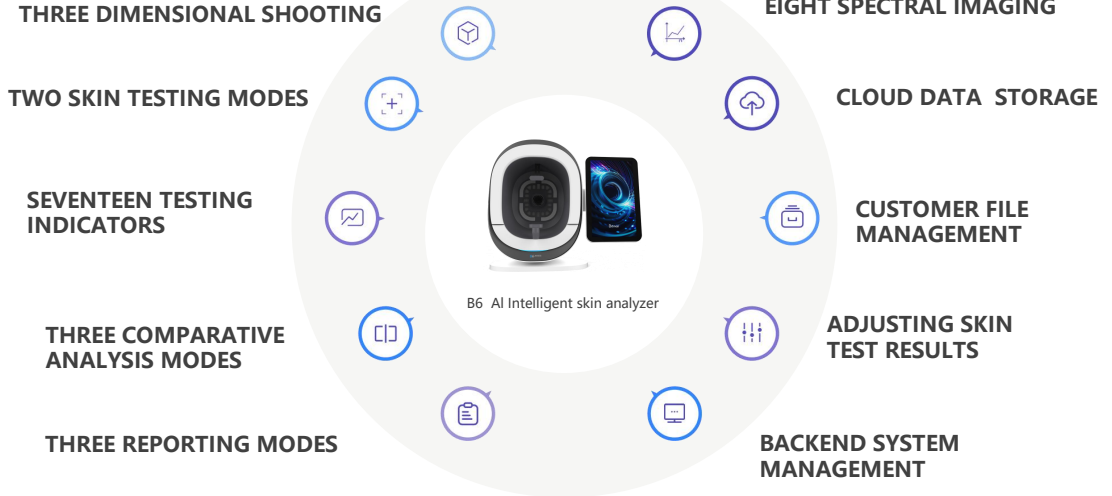

Three dimensional shooting

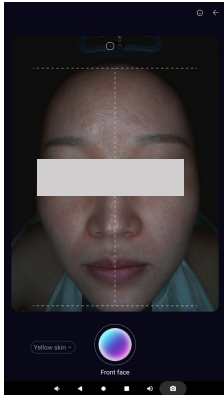

Front face

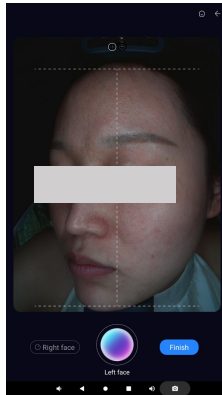


Left face



Right face

Four skin tones available

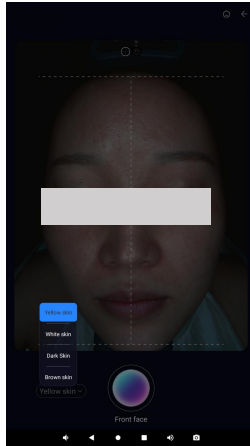

Yellow Skin

White skin

Dark skin

Brown skin


- Eight spectral images




Seventeen testing indicators



Three comparison modes

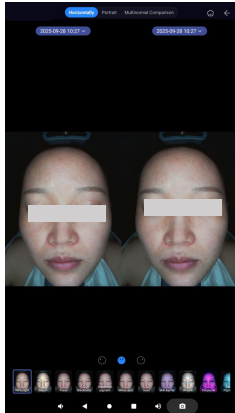


Parallel comparison

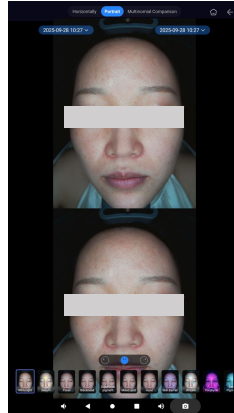

Vertical comparison


Three comparison modes

Comparison of the effects of single indicators before and after.

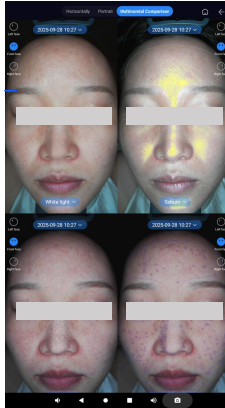


Comparison of four figures

Analysis and comparison of multiple problem indicators before and after skin care.

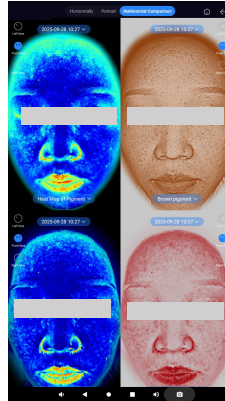

Comparison of four figures

Three reporting modes-single independent report 1

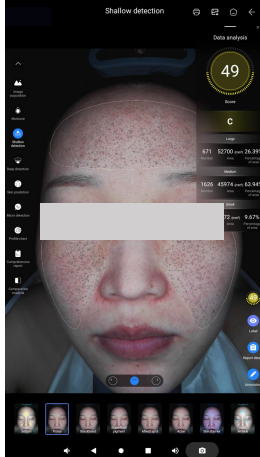

Introduction to data analysis

Data score



Multidimensional perspective

Multi-angle comprehensive image display


- 1.Score
- 2.Level
- 3.Number
- 4.Area
- 5.Percentage of area

Let consumers accurately understand their skin problems

Accurately quantify the underlying effect and empower doctors to treat



It is divided into 5 levels according to the skin condition from high to low and marked with different colours.

A Green B Blue C Yellow D Orange E Red

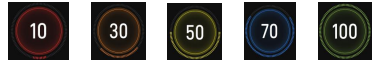

E D C B A

Three reporting modes-single independent report 1

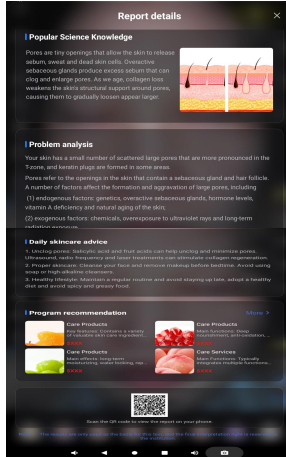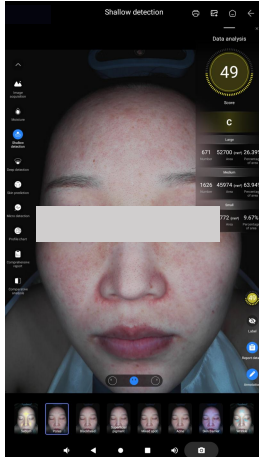


Popular Science Knowledge

Problem analysis

Daily skincare advice

Program recommendation

The classification of 16 skin types in Baumann Skin Typing System

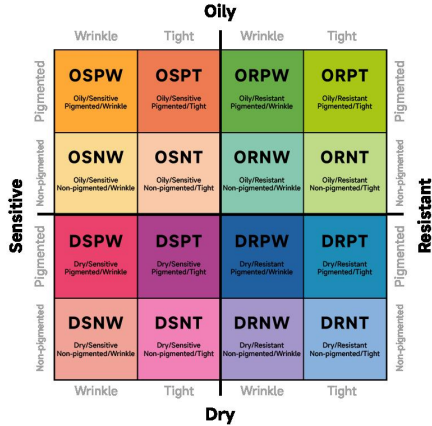

Type	English	Type	English
DSPT	Dry, Sensitive, Pigmented, Tight	OSPT	Oily, Sensitive, Pigmented, Tight
DSPW	Dry, Sensitive, Pigmented, Wrinkle	OSPW	Oily, Sensitive, Pigmented, Wrinkle
DSNT	Dry, Sensitive, Non - Pigmented, Tight	OSNT	Oily, Sensitive, Non - Pigmented, Tight
DSNW	Dry, Sensitive, Non - Pigmented, Wrinkle	OSNW	Oily, Sensitive, Non - Pigmented, Wrinkle
DRPT	Dry, Resistant, Pigmented, Tight	ORPT	Oily, Resistant, Pigmented, Tight
DRPW	Dry, Resistant, Pigmented, Wrinkle	ORPW	Oily, Resistant, Pigmented, Wrinkle
DRNT	Dry, Resistant, Non - Pigmented, Tight	ORNT	Oily, Resistant, Non - Pigmented, Tight
DRNW	Dry, Resistant, Non - Pigmented, Wrinkle	ORNW	Oily, Resistant, Non - Pigmented, Wrinkle

Three reporting modes-Comprehensive analysis report 2



Personal information

16-type skin classification and comprehensive score


On the left is the skin water content and on the right is a single report.

Overall recommendation

Single indicators below C will be realised in the form of radar charts.

Real-time analysis of problems and nursing suggestions based on Bauman's 16 types of skin classification

Three reporting modes-Comprehensive analysis report 2

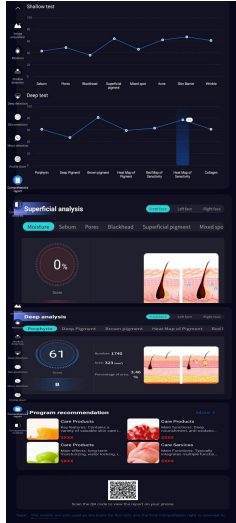

Shallow test


Comparative analysis,
only compare and
analyse with the results
of the last shooting.

Deep test

Shallow analysis

Deep analysis


Download the report

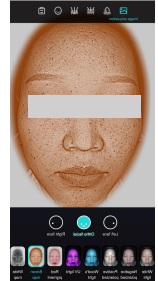

H5 Mobile phone report

Scan the code on your mobile phone to get the report

Three kinds of reporting modes Innovation Independent Editing Report 3

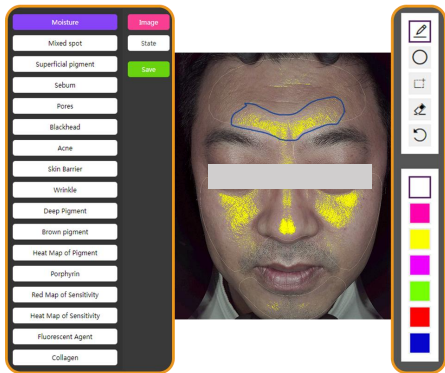


The background automatically selects the image to be edited for annotation

Customize any test results you want

You can debug results for each indicator

Free copy editing

Parameter adjustment

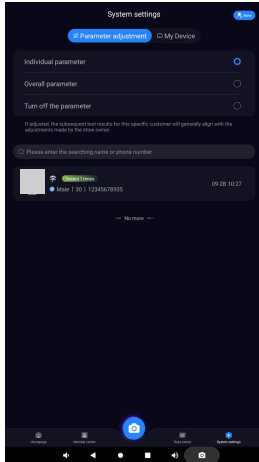


Parameter adjustment




Data result optimization manual debugging

It can be adjusted as a whole

My device

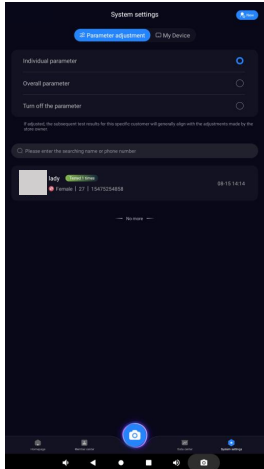


The language
can be adjusted

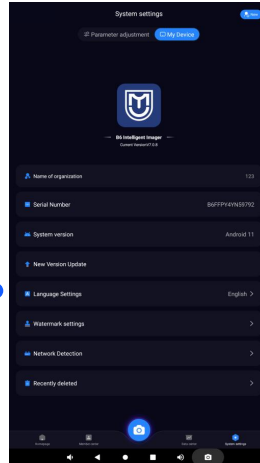


Serial number

You can find the
instrument
problem through
the background
and solve it.


Image export function

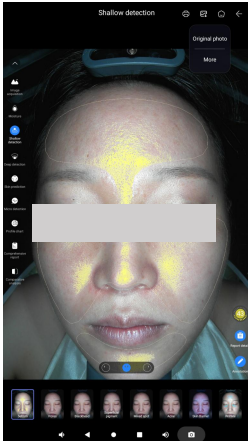

Exporting images

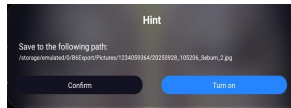


This figure is exported.

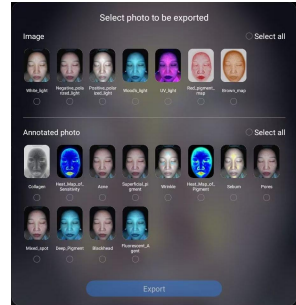

Export multiple images


Member Center

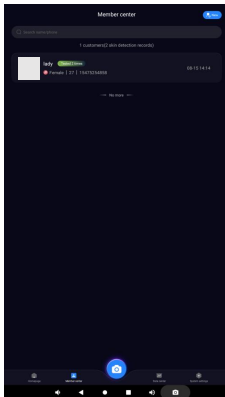

File Management

The test items with scores below 60 will be listed here.

Bauman 16-type classification

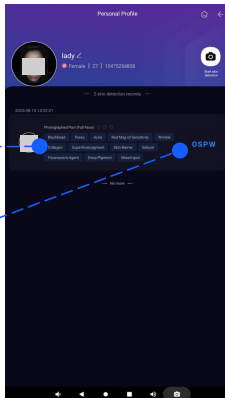

Profile

Cloud storage

One click search

Number of skin tests

Data Center

Choose the year and the month

The number of people entering the store in a month or a year can be displayed more intuitively through the line chart

Ratio of men to women

He men and women who enter the store can see the ratio of male and female women who enter the store through the data form more clearly and intuitively

Problem proportion

In the form of data, the current population is divided into five levels with different colors to represent the data in each indicator

Available year

Age distribution

Through the age division of the number of people entering the store in the form of data.



Backend management



**Recording
nursing
project freely**

**Unified
management
of multiple
equipments**

**Customer
profile
management
in real-time**

**Review
and edit
detection
record in
real-time**

Lead technology create beauty



B6 AI Intelligent skin analyzer

02

PARAMETERS


Parameters

Parameter information

B6

Type number	B6
Image element	Industrial grade 36 million
Light spectrum	White light, positive polarization, negative polarization, Wu, UV, erythroid, brown, white
Shading mode	Pull-down hood
Product material	Industrial grade ABS
Electric source	AC100-260V,50/60Hz
Product size	695 x 305 x 569mm (after folding: 516x 402 x 569mm)
Central processing unit	The RK3568 quad-core Cortex-A55 has a maximum frequency of 2.0GHz
Main plate	R10-S6810 motherboard with Android11 system
Internal memory	Dual channel LP DD R4, 4G
Hard disk	MMC5.1, 32G
Mode of operation	Multipoint capacitive touch
Screen feature	It folds 90 degrees
WIFI	Built-in dual-band WIFI(2.4G, 5G)
HDMI	1 Pcs
USB	2 Pcs
Screen scale	16:9
Screen size	15.6 inches
Screen pixel	1920 * 1080

Hardware parameters

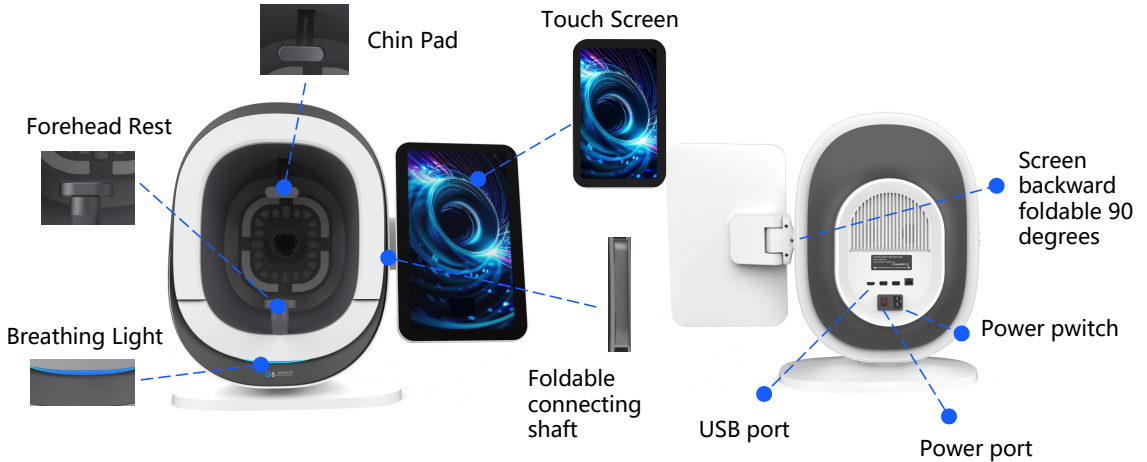

【B6 AI Intelligent skin analyzer】

Lead technology create beauty






B6 AI Intelligent skin analyzer

03

EIGHT
SPECTRAL
IMAGE
ANALYSIS


Microscopic detection

- Microscopic imaging of skin and scalp details issues

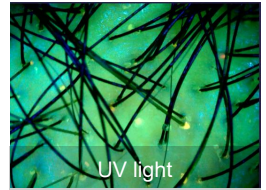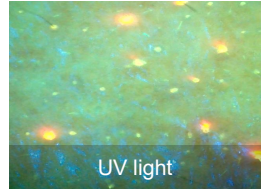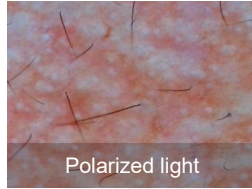

Multiple light source targeting analysis

Eight spectral image analysis

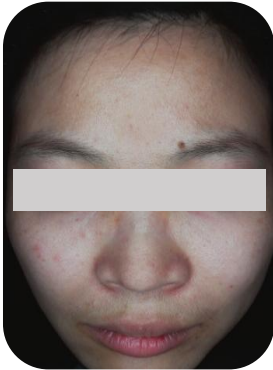

White light

THEORY

Visible spots and other blemishes on the skin surface (acne, spots, wrinkles, pores, etc.) under natural light sources, which are mainly used as the basis for other spectral image comparison.



Acne



Dark spots




Wrinkle



Pore

Eight spectral image analysis

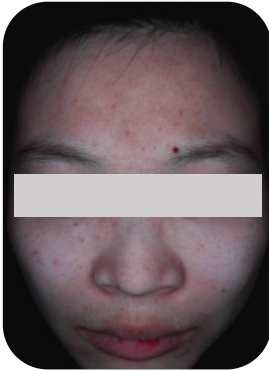


Negative polarized light

THEORY

Using negative polarized technology to filter out the refracted light on the skin surface, so that you can clearly examine the light brown, tan, dark brown, light yellow or dark red skin lesions; It can distinguish the condition of capillaries, facial acne, uniformity skin and other skin problems.

Dark spots

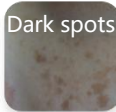

Pigment



Acne




Eight spectral image analysis




Positive polarized light

THEORY

Positive polarized light can improve the clarity of superficial texture, magnify local details, so as to clearly observe the smoothness of skin, fine lines and wrinkles and bumps (wrinkles, pores, Acne scars, Acne, etc).

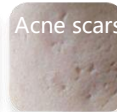

Eight spectral image analysis

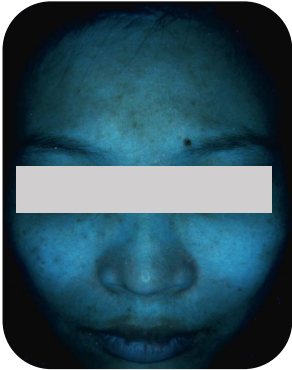

Wood' s light

THEORY


Wood' s light can detect deep pigments in dermis. The principle behind this is that melanin does not fluoresce after exposure to ultraviolet radiation, allowing melanin to stand out more clearly with stronger contrast.

Spots

Spots

Fluorescence

Eight spectral image analysis

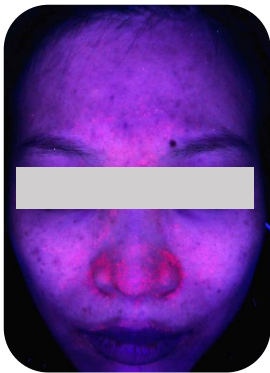


UV light

THEORY

Under UV light source, the content and distribution of the purple pigment bilirubin are displayed clearly through fluorescence, which can be used for the auxiliary diagnosis and efficacy observation of pigmentary dermatoses, pore issues, skin infections, and porphyria.



Eight spectral image analysis

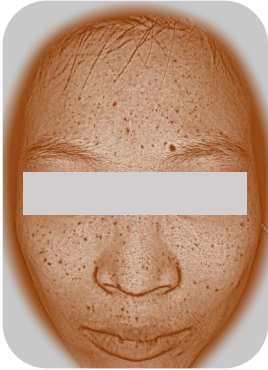


Brown light

THEORY

The position, area, shape, and severity of subcutaneous facial UV spots are processed by using RBX light source technology, which demonstrate skin damage from UV radiation and the accumulation of subcutaneous melanin.



Eight spectral image analysis

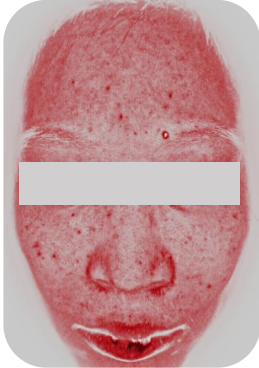


Red light

THEORY

Using the RBX light source technology, it is applied for analyzing the subcutaneous hemoglobin on the face as well as inflammatory pigmentation (such as sensitive skin, lesions, acne, red blood vessels, etc., which are skin inflammatory conditions)

Acne

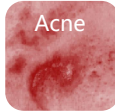

Blood
streak



Mole



Eight spectral image analysis



Mixed light

THEORY

Skin texture roughness and collagen loss were revealed by polarizing analysis.

Rough
texture

Wrinkle




B6 AI Intelligent skin analyzer

04

24

**DETECTION
FUNCTION**

24 Detection function- Annotation function 1

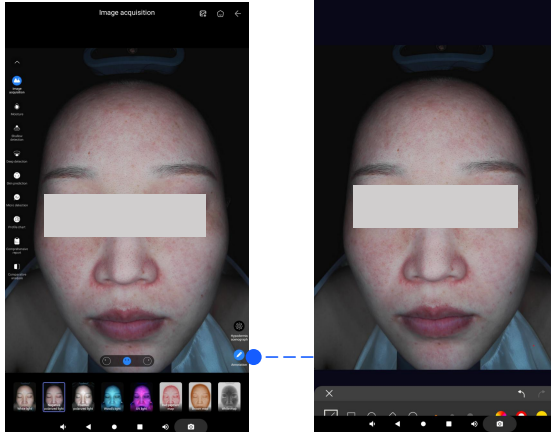

THEORY

Annotation function uses tools to mark visible skin issues (e.g., redness, spots) on images, intuitively showing conditions and facilitating skincare advice explanation.

24 Detection function-Moisture test report 2



Zoned moisture detection

Partition detection of skin and facial moisture

Skin water content score




Moisture detection handle

The water content of the skin is sorted from high to low according to five levels and marked with color.



23 Detection function-Moisture test report 2

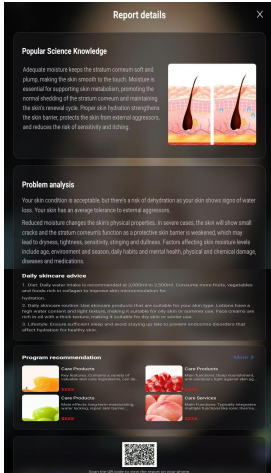


Popular Science Knowledge

Problem analysis


Daily skincare advice

Program recommendation


Scan the code to view the report

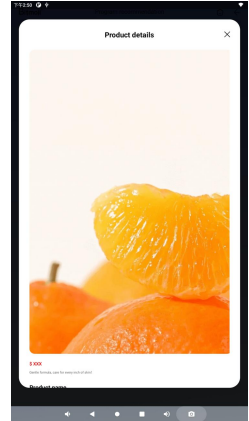

24 Detection function-Under - skin perspective

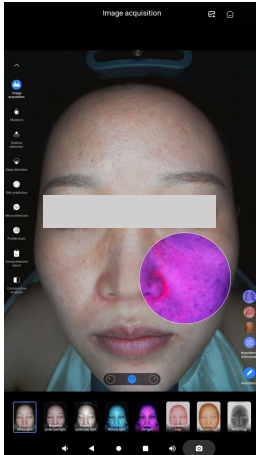


Under - skin
perspective

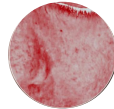

Clicking the "Muscle Layer Perspective" button on the right side will reveal three icons: the UV light map, the melanin map, and the brown map. You can freely move the perspective lens to compare issues such as pigmentation, red blood vessels, and skin inflammation on the superficial layer and the base layer of the skin.

24 Detection function

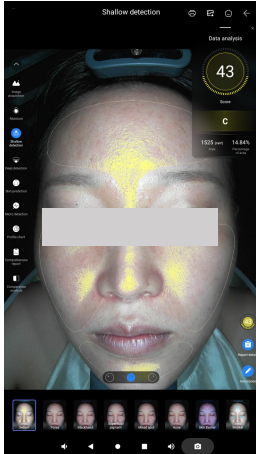


Sebum

Image Analysis

- ❑ The oil secretion of the skin surface can be checked under positive polarized light source.
- ❑ The algorithm displays areas of the skin with active oil secretion through yellow ,through the form of data, you can see the oiliness of facial skin more clearly and intuitively.
- ❑ Excess oil is one of the factors that trigger acne growth, so please take good oil control care if you have acne.

24 Detection function

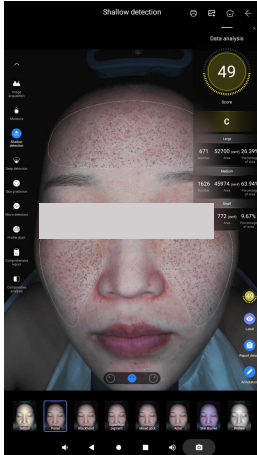

Pores

Image Analysis




- ❑ Under negative polarized light source, it is possible to check if enlarged pores have formed on the skin surface.
- ❑ The algorithm uses the RBX technology to process and present the areas with large pores in the skin as dark gray small circles; it can also show the situation of already enlarged pores in the facial skin in a clearer and more intuitive data form.
- ❑ Pore clogging refers to the pores on the surface of the skin being blocked, which prevents sebum from being discharged normally, accompanied by the accumulation of stratum corneum and dirt. this phenomenon usually manifests itself in the form of blackheads, whiteheads or acne, and in severe cases may lead to skin problems such as acne and folliculitis.

[illegible]

Image Analysis

- ❑ Under negative polarized light source, you can see blackheads formed by pores clogged by oil in the T-zone.
- ❑ The algorithm uses the RBX technology to process and present the area of blackheads in the T-zone as small dark gray circles; it enables a clearer and more intuitive view of the blackheads on the nose area through the data.
- ❑ Blackheads are formed by excess oil accumulation in the nose area of the skin and air oxidation. Areas with large pores are more likely to accumulate and store oil and dust in the air, so it is necessary to clean and moisturize in time to reduce the formation of large pores.

24 Detection function

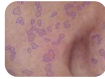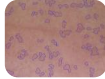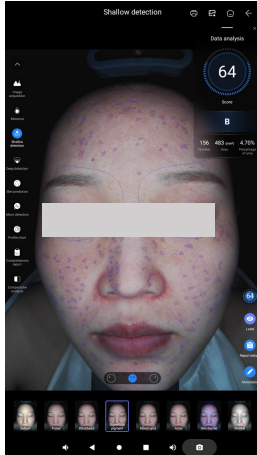


Superficial Pigment

Image Analysis

- ❑ Superficial pigment refers to pigmentation that has formed on the superficial layer of the skin, including: acne scars, spots, inflammatory pigmentation, etc.
- ❑ The coverage of pigmentation may exist in both deep and shallow layers. you can compare the image with the deep pigment image. If the shallow layer shows pigmentation but the deep layer shows no pigmentation, it means that the pigment is only deposited in the superficial layer of the skin.
- ❑ The algorithm marks the pigmented area with a purple polygon curve, and the shallow pigment can be seen more clearly and intuitively through the form of data.

24 Detection function

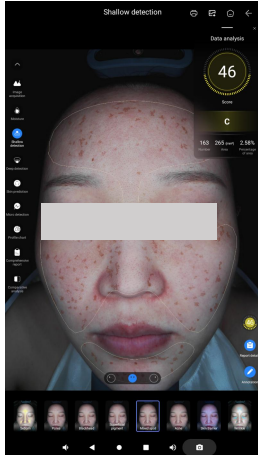


Mixed Spot

Image Analysis

- ❑ Under negative polarized light, we can see the distribution of mixed spots on the skin surface. The algorithm identifies the facial complex spot area and marks it with a brown block.
- ❑ The mixed spot map shows skin pigmentation such as melasma, age spots, and freckles. Melasma is a darker patch on the skin that can appear brown, black, or dark brown. Melasma may expand over time, especially if daily sun protection and skin care are not taken care of. Some melasma may be slightly raised and feel slightly convex to the touch.
- ❑ The algorithm marks the mixed spot area with brown color blocks, and the mixed spot situation can be seen more clearly and intuitively through the form of data.

24 Detection function

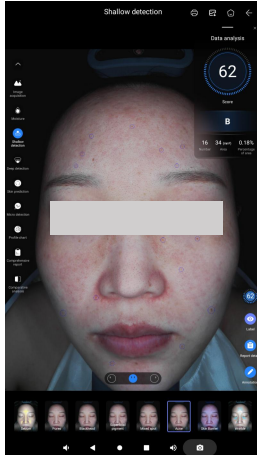


Acne

Image Analysis

- ❑ Look at the distribution of skin acne and superficial redness under negative polarized light.
- ❑ When pores are clogged with oil and dust, it is easy to fester or form inflammation, which will then turn into acne and acne.
- ❑ The algorithm identifies the distribution area of facial acne and marks it with blue circles. The more the number and the more obvious the redness of the skin, the more serious the skin acne problem is, and the skin needs to be oil-controlled to unclog the pores and eliminate inflammation. You can see the acne situation more clearly and intuitively through the form of data.

24 Detection function



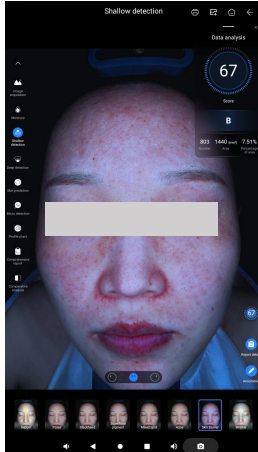


Skin Barrier

Image Analysis

- ❑ We can check the skin barrier health under negative polarized light source.
- ❑ The barrier image shows the skin redness problem and the distribution of red blood streaks. The formation of red blood streaks is mainly due to the damage of keratin, the weakness of the epidermis, and the long-term damage of the capillary position, which leads to vascular dilation and congestion.
- ❑ The redder the area, the more severe the damage to the skin barrier. This can be used as a reference for judging the sensitivity of the skin and the presence of inflammatory areas. The damaged condition of the barrier can be seen more clearly and intuitively in the form of data.

24 Detection function

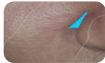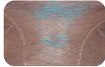
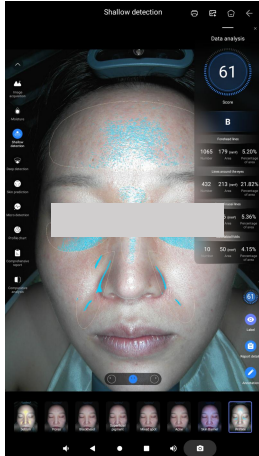


Wrinkle

Image Analysis

- ❑ The texture of the skin surface can be viewed under a positive polarized light source.
- ❑ The wrinkle image shows the roughness of the skin texture, such as large pores, dry lines, fine lines, and wrinkles. It can be used as a reference for judging the fineness of the skin and the loss of collagen.
- ❑ The algorithm identifies the patterns on the facial skin and marks the distribution of skin wrinkles in five areas (forehead wrinkles, nasal bridge wrinkles, eye area wrinkles, crow's feet wrinkles, and laugh lines) with Lake blue dotted lines. The more discontinuous the lines, the rougher the skin. This way, the wrinkles can be seen more clearly and intuitively in a data form.

24 Detection function


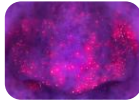
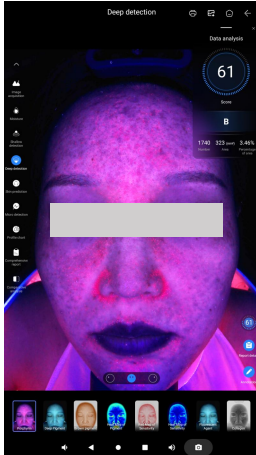


Porphyrin

Image Analysis

- ❑ The brick-red fluorescent spots in the picture are propionibacterium acnes and malassezia. these two bacteria will aggravate the occurrence of skin acne, so they can be used as a basis for judging skin acne. Through the form of data, the situation of porin can be seen more clearly and intuitively.
- ❑ The living environment of propionibacterium acnes and malassezia must have oil, so they can be used as a basis for judging the accumulation of oil in skin pores.

24 Detection function

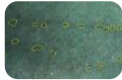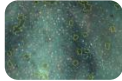
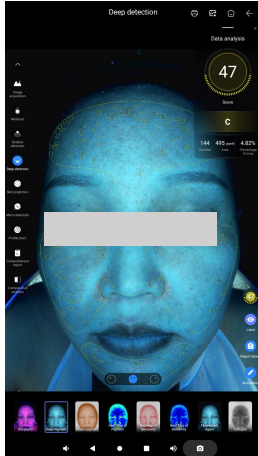


Deep Pigment

Image Analysis

- ❑ The bottom part of the figure shows wood's light. The yellow area represents the comprehensive color spot area identified by the algorithm and is marked with polygonal curves. This way, the situation of deep pigmentation can be seen more clearly and intuitively in the form of data.
- ❑ The dark (black, brown) block or dot skin that appears on the face is a display of skin pigmentation (such as melasma, freckles, malar spots, inflammatory pigmentation, acne marks, hemoglobin aggregation, etc.).
- ❑ The pigmentation in the deep layer of the skin can be compared with the sensitivity to determine whether it is an inflammatory hemoglobin accumulation or a spot problem.

24 Detection function



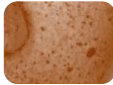
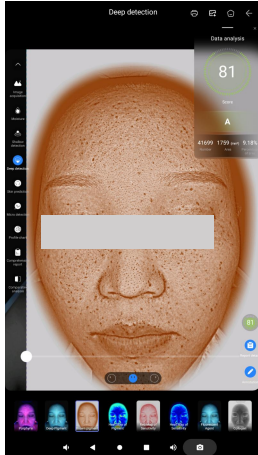


Brown Pigment

Image Analysis

- ❑ The depth of the overall brown color of the skin is mainly related to the skin color. People with darker skin or more hemoglobin have darker overall pigmentation.
- ❑ The areas with heavier pigmentation in the image are mostly those with higher pigment concentration density.
- ❑ Through the form of data, the brown pigment can be seen more clearly and intuitively.

24 Detection function

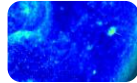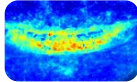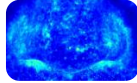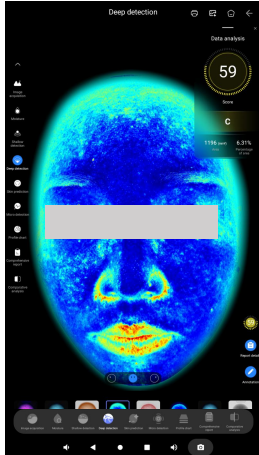

Heat Map of Pigment


As shown in these images

- ❑ Heat map of pigment can check the distribution of pigment deep in the skin.
- ❑ The algorithm identifies the distribution of pigments on the face and presents it in the form of a heat map. Different colors are used to represent the distribution of spots, moles, and scars visible to the naked eye under negative polarized light. Red indicates severe skin pigmentation, yellow for medium, green for lighter skin, and blue for normal skin. The situation of pigmentation can be seen more clearly and intuitively through the form of data.
- ❑ Pigment production mechanism: the body's own regulation, physical or chemical factors stimulate melanocytes, so that their number increases and activity increases. The melanin produced cannot be completely removed with the stratum corneum and blood circulation, and finally deposits in the local skin.

- The algorithm identifies the distribution of pigments on the face and presents it in the form of a heat map. Different colors are used to represent the distribution of spots, moles, and scars visible to the naked eye under negative polarized light. Red indicates severe skin pigmentation, yellow for medium, green for lighter skin, and blue for normal skin. The situation of pigmentation can be seen more clearly and intuitively through the form of data.

- ❑ **Pigment production mechanism:** the body's own regulation, physical or chemical factors stimulate melanocytes, so that their number increases and activity increases. The melanin produced cannot be completely removed with the stratum corneum and blood circulation, and finally deposits in the local skin.




24 Detection function

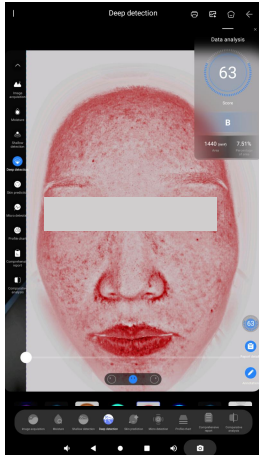





Red Map of Sensitivity

As shown in these images

- ❑ In the Red Map of Sensitivity, we can observe the redness of the skin's surface layer and the distribution of red blood vessels.
- ❑ The picture clearly shows the distribution of red blood vessels, indicating that the skin is thin and sensitive, and requires proper protection and care.
- ❑ The shade of the background color in the red photo is related to the overall skin tone. Those with less hemoglobin will have a lighter color.
- ❑ The areas with a higher concentration of red color indicate a more concentrated accumulation of hemoglobin in the skin, which can be used as a reference for judging the sensitivity of the skin and the presence of inflammatory areas. The degree of skin sensitivity can be clearly and intuitively seen through the data form.

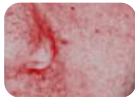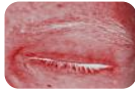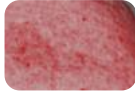

24 Detection function

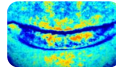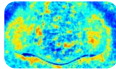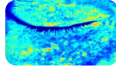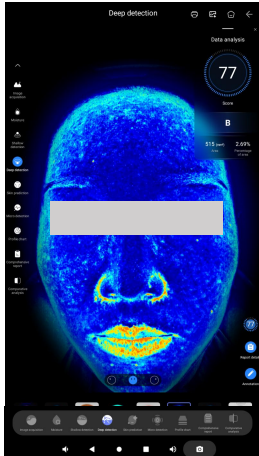


Heat Map of Sensitivity

As shown in these images

- ❑ The heat map of sensitivity represents skin sensitivity. When the skin shows significant redness and thinning of the stratum corneum, it becomes more susceptible to external stimuli and damage, leading to issues such as dryness, sensitivity, and redness.
- ❑ The heat map of sensitivity is based on the distribution of subcutaneous capillaries, with areas of greater sensitivity having more capillaries. Visible redness and acne on negative polarized light images indicate areas of severe sensitivity.
- ❑ The algorithm uses different colors to indicate varying degrees of sensitivity and their distribution on the skin. Areas with severe sensitivity are shown in deep red, including the lips; medium sensitivity is represented in yellow, mild sensitivity in green, and normal skin appears in blue. The sensitivity is more clearly and intuitively reflected in the form of data.

24 Detection function



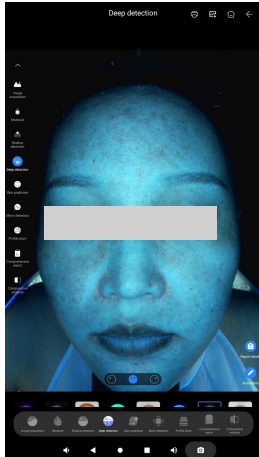


Fluorescent Agent

As shown in these images

- ❑ Fluorescent agent and pigments might both appear in facial imaging. To assess the fluorescent agent, focus specifically on the fluorescence response.
- ❑ The difference between fluorescent agent and porphyrins is as follows porphyrins exhibit brick-red fluorescent agent, while fluorescent agent display intense blue light and usually appear as large, sheet-like areas.
- ❑ The difference between fluorescent agent and facial .
- ❑ Dust is as follows facial dust appears as white, bright, floating, and short, wispy lines on the surface, while fluorescent agent typically display bright colors and are often more diffuse or spread over larger areas.

24 Detection function

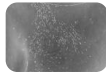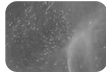
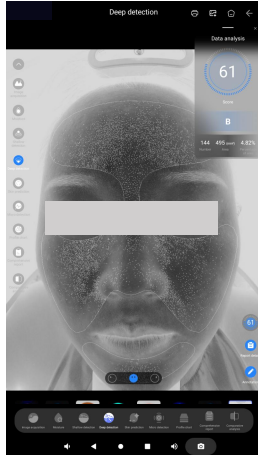


Collagen

As shown in these images

- ❑ Under the mixed light image, we can observe the loss of collagen on the skin surface.
- ❑ The mixed light image shows a situation where the skin texture is rough. For example: large pores on the skin, fine lines and wrinkles. It can be used as a reference for judging the skin's smoothness and the loss of collagen.
- ❑ The white dotted lines in the mixed light image indicate the loss of skin collagen and the rupture of elastic fibers. This way, the loss of collagen can be clearly and intuitively observed in a data form.




24 Detection function

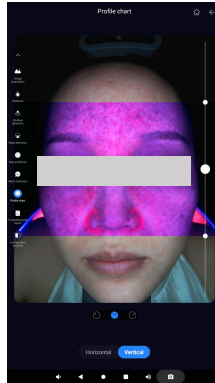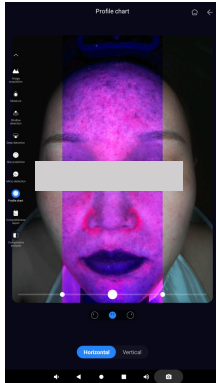

Profile chart

Through white light, negative polarized light and UV light source comparison, multi-dimensional, deep analysis of skin problems.



24 Detection function

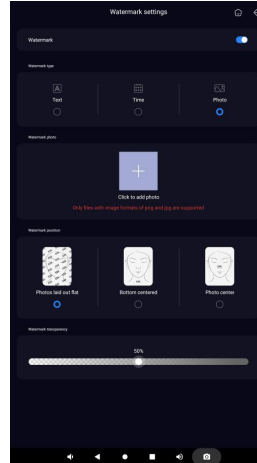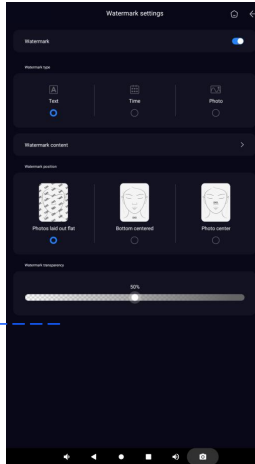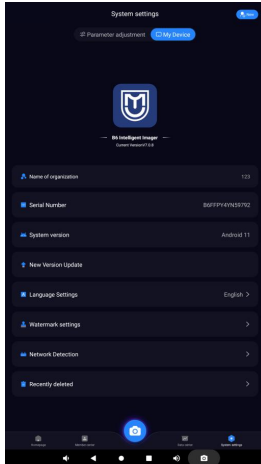

Watermark function

Customizable text
watermark, image
watermark


24 Detection function

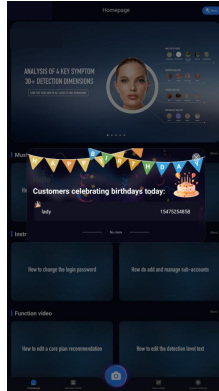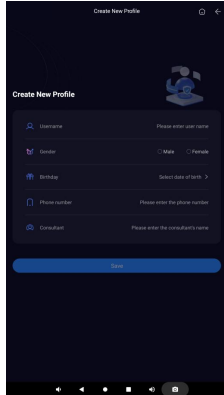

Birthday reminder function

Automatically send out the list of customers whose birthdays are on that day, in order to enhance customer retention.

24 Detection function

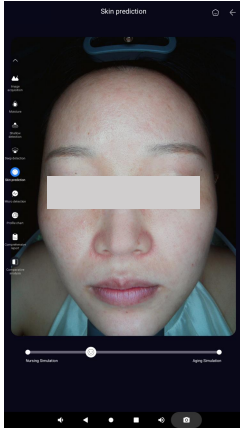


Nursing simulation

Skin prediction

Deeply predict the future of the skin and awaken customers' desire for young skin.

According to the current skin condition of the customer, through the training of AI large model, simulate the skin condition after customer care and the aging situation of different age groups.

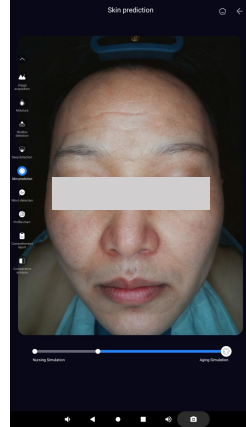

Aging simulation

【B6 AI Intelligent skin analyzer】

Lead technology create beauty




B6 AI Intelligent skin analyzer

05


**MICROSCOPIC
DETECTION
SKIN DETAILS**

Microscopic detection - Skin Details

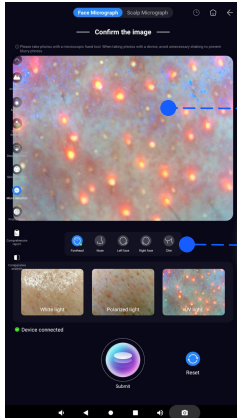

Local microscopic display

Location selection, partition detection.

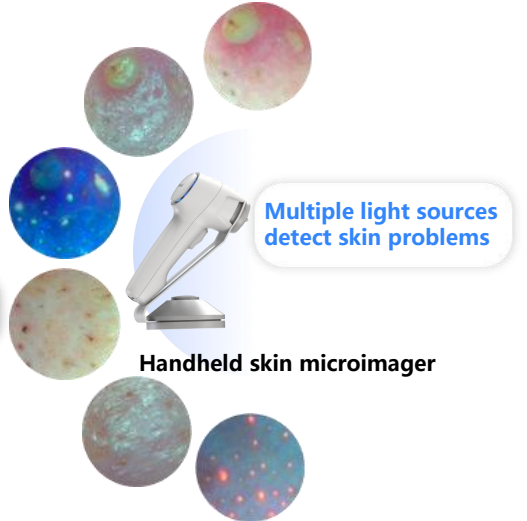


Multiple light sources detect skin problems

Handheld skin microimager

Microscopic detection - Skin Details




Independent
record tracking
for each
customer

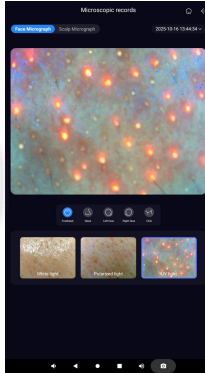

Comparison
display for the
same area

Microscopic examination - scalp details

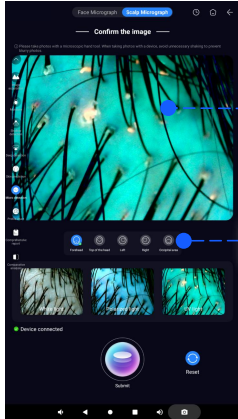

Local microscopic display

Location selection, partition detection.




Multiple light sources
Discover scalp problems

Handheld skin microimager



【B6 AI Intelligent skin analyzer】

Lead technology create beauty

Microscopic examination - scalp details

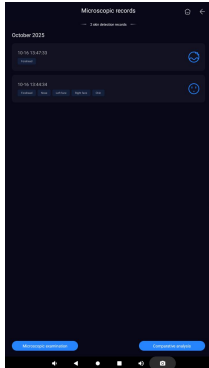

Independent
record tracking
for each
customer

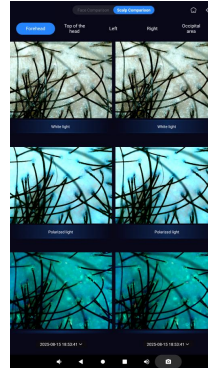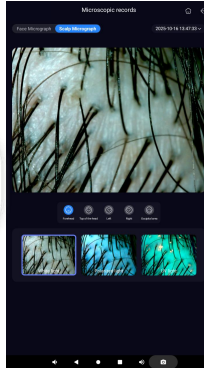

Comparison
display for the
same area

Welcome you join us!